

# INSTALLATION, OPERATION AND MAINTENANCE MANUAL

CLIMATE CONTROL DEVICE 4000-01



# TYPE CE, CONTINUOUS FLOW, SUPPLIED-AIR RESPIRATOR THESE RESPIRATORS ARE APPROVED ONLY IN THE FOLLOWING CONFIGURATIONS

|         |                 |                       |                      |        |  |         |                    |         | RE                           | ESI     | PIR      | RA                 | ГО     | R C    | CO   | MF     | 201              | NE     | NΊ     | ۲S     |                               |        |        |        |        |        |        |        |        |            |
|---------|-----------------|-----------------------|----------------------|--------|--|---------|--------------------|---------|------------------------------|---------|----------|--------------------|--------|--------|--|--------|------------------|--------|--------|--------|-------------------------------|--------|--------|--------|--------|--------|--------|--------|--------|------------|
| TC No.  | Protection<br>1 | M<br>O<br>D<br>E<br>L | Alternate<br>Helmets |        | Alternate<br>Flow<br>Regulator<br>and Belt |         | Alternate<br>Capes |         | Alternate Quick<br>Disconect |         |          | Alternate<br>Hoses |        |        | Breat-<br>hing<br>Air<br>Tube Alternate<br>Visor |        | Alternate Lenses |        |        |        | Cautions and<br>Limitations 2 |        |        |        |        |        |        |        |        |            |
|         |                 |                       | NV2001               | 02-803 | NV2016                                     | 4000-01 | NV2015             | 4000-20 | NV2002                       | NV2002L | NV2002XL | 02-813             | NV2025 | NV2024 | NV2032   | NV2034 | NV2035           | NV2028 | NV2029 | NV2035 | NV2036                        | NV2021 | NV2004 | 02-804 | NV2018 | NV2031 | NV2017 | 02-810 | 02-811 |            |
| 19C-363 | SA/CF/SB        | NOVA<br>2000          | Х                    |        | Χ  | Χ       |                    | Χ       | Χ                            | Х       | Χ        |                    | Χ      | Х      | Χ  |        |                  | Х      | Х      |        |                               | Х      | Х      |        | Х      | Х      | Х      |        |        | ABCDEJMNOS |
| 19C-417 | SA/CF/SB        | ASTR0                 |                      | X      | Χ  |         |                    |         |                              |         |          | Χ                  | Χ      | Χ      | Χ  |        |                  | Х      | Х      |        |                               | Х      |        | X      |        |        |        | X      | Х      | ABCDEJMNOS |
| 19C-422 | SA/CF/SB        | NOVA<br>2000          | X                    |        |  |         | Χ                  |         | Χ                            | X       | Χ        |                    |        |        |  | Χ      | Χ                |        |        | Χ      | X                             | X      | Х      |        | X      | Х      | Х      |        |        | ABCDEJMNOS |
| 19C-423 | SA/CF/SB        | ASTR0                 |                      | X      |  |         | Χ                  |         |                              |         |          | Χ                  |        |        |  | Х      | X                |        |        | Χ      | X                             | Х      |        | Χ      |        |        |        | X      | X      | ABCDEJMNOS |

#### 1. PROTECTION

CF- Continous flow

SA- Supplied Air

SB- Supplied - Air Abrasive Blast

#### 2. CAUTIONS AND LIMITATIONS

- A- Not for use in atmospheres containing less than 19.5 percent oxygen.
- **B-** Not for use in atmospheres immediately dangerous to life and health.
- **C-** Do not exceed maximum use concentrations established by regulatory standards.
- **D-** Air-line repirators can be used only when the repirators are supplied with respirable air meeting the requirements of CGA G-7.1 Grade D or Higher Quality.
- **E-** Use only the pressure ranges and hose lenghts specified in the User's instructions.
- **J-** Failure to properly use and maintain this products could result in injury or death.
- M- All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N- Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- **0-** Refer to User's instructions, and/or maintenance manuals for information on use and maintenance of these repirators.
- **S-** Special or critical User's instrucions and/or specific Limitations apply. Refer to User's instrucions.

# INTRODUCTION

The 4000-01 Cold Air Tube is NIOSH approved to supply a continuous flow of cool compressed air to the Supplied Air Respirator. All Genuine Parts must be present and correctly assembled to constitute a NIOSH approved respirator.

#### **AIR QUALITY**

The respirator must be supplied at all times with clean breathable air Grade D or better. Always use an inlet filter on your compressor - locate your compressor in a clean environment. The AIR QUALITY must be constantly monitored to ensure Grade D air is supplied at all times. A carbon monoxide alarm must always be used.

# **⚠ WARNING ⚠**

Do not connect the air supply hose to nitrogen, toxic gases, inert gases or any other non Grade D air source.

#### **AIR PRESSURE**

The air pressure supplied to the respirator must be continuously monitored to ensure it keeps within the specified range as detailed on the breathing air pressure tables.

# **△** WARNING **△**

Make sure the correct table is used. When using 4000-01 use (Fig 3)

Breathing Air Pressure Table.

# **△** WARNINGS **△**

- **DO NOT** use the 4000-01 until you have been trained in its use, maintenance and limitations by a qualified individual (appointed by your employer) who has extensive knowledge of the 4000-01
- **DO NOT** ingest the ice that may form around the cold air outlet as it could contain oils or toxic element that could cause illness or death.
- USE ONLY NIOSH approved NOVA 2000 air supply hose.
- DO NOT modify or alter parts of the respirator. The use of non approved parts voids the entire approval of the Respirator.

# SET UP AND USE

- 1. To supply cold air to the respirator connect the outlet to the NV2021 breathing tube and tighten by hand. (Fig 1).
- 2. Connect the NIOSH approved air supply hose to a breathing air source, supplying Grade D orbetter quality air. Connect the air supply hose to the 4000-01 quick disconnect tail. (Fig 2).
- With air flowing through the respirator, adjust the air pressure at the point of attachment. (Fig 2) to the recommended air pressure specified in the BREATHING AIR PRESSURE TABLE (Fig 3).


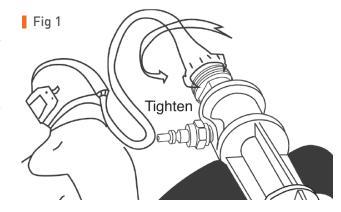


Fig 2

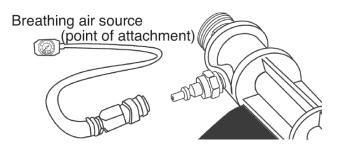

# \Lambda WARNING \Lambda

Make sure you understand the breathing air pressure table before using the respirator. Make sure you set the correct air pressure for the correct FIG 2 supply hose length.

Fig 3

|                                      |                            | BREATHING AIR I   | PRESSURE TABLE                              |                                 |   |
|--------------------------------------|----------------------------|-------------------|---|---------------------------------|---|
| AIR<br>SOURCE                        | BREATHING TUBE<br>ASSEMBLY | AIR SUPPLY HOSE   | AIR SUPPLY HOSE<br>LENGTH                   | MAX NUMBER OF<br>SECTIONS       | PRESSURE RANGE<br>(PSIG Air)  |
| Portable or Stationary<br>Compressor | 4000-01<br>Cold Air Tube   | NV2028/<br>NV2029 | 25<br>50<br>100<br>150<br>200<br>250<br>300 | 1<br>1<br>2<br>3<br>4<br>5<br>6 | 55 - 56<br>55 - 57<br>60 - 62<br>65 - 67<br>70 - 72<br>77 - 78<br>80 - 82 |

4. To obtain cooler air turn the regulator control knob anti cockwise, this will increase the airflow out of the exhaust port (Fig 4). Turn the knob clockwise to obtain air closer to ambient temperature.


## 

DO NOT USE THE 4000-01 when ambient temperature is below 68; (20;C) as ice could form in the cold air outlet and reduce airflow.


5. With the compressed air flowing into the respirator at the correct pressure range for the correct supply hose length you are now ready to enter the work area. Do no remove your respirator until you are outside the blasting area.



Leave work area immediately if you are unwell, airflow stops or slows down, breathing becomes difficult, you become dizzy, too cold or if any respirator part is damaged.

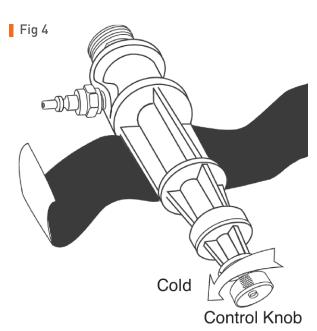

# **CCT 4000 PARTS LIST**

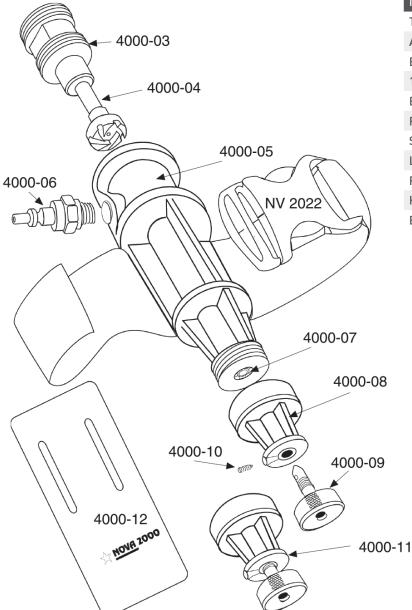

| Item                       | PART No. |
|----------------------------|----------|
| Top Threaded Cap           | 4000-03  |
| Air Spinner                | 4000-04  |
| Body                       | 4000-05  |
| 1/4" Quick Disconnect Tail | 4000-06  |
| Breaker                    | 4000-07  |
| Regulator Body             | 4000-08  |
| Spindle & Knob             | 4000-09  |
| Locking Screw              | 4000-10  |
| Regulator Assembly         | 4000-11  |
| Heat Sheild                | 4000-12  |
| Belt                       | NV2022   |

# **△ WARNING △**

Due to the high noise levels created during abrasive blasting earplugs must be worn at all time. During Abrasive Blasting operations do not use Silica Sand, the use of Silica Sand can cause serious respiratory disease or death, refer: www.osha.gov.Silica.IT69D\_1.html

| NOTES |  |
|-------|--|
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |
|       |  |

| NOTES |
|-------|
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |
|       |

# OTHER PRODUCTS









ISO9001
ACCREDITED COMPANY

# **DID YOU KNOW**

We can also supply you with the following...

#### **NOVA 2000 BLASTING HELMET**

The NOVA 2000 has been designed for blasting professionals.

The NOVA 2000 is the most advanced abrasive blasting helmet available, with fully padded pillow foam liner for extra comfort and noise attenuation, wide angle vision, even air distribution assures a pleasant working environment, optional climate control devices. Things cannot get much better.

### **ASTRO BLASTING HELMET**

ASTRO is your number one choice economic Abrasive Blasting Helmet which meets safety standards worldwide.

This is the lightweight, robust helmet you have been waiting for. It helps to reduce fatigue, reduce down time and increase productivity.

#### **HEAD SOCKS & SPRAY HOODS**

Stay cleaner with Head Socks and Spray Hoods. One size fits all. The cool comfortable Spray Hood protects from over spray and sanding dust. The Head Sock can be worn as a hygienic liner.

#### **HEAD SOCKS & SPRAY HOODS**

The RADEX Airline Filter offers you more capacity, versatility & filtration. You can combine the versatility of either floor or wall mounting with increased filtration capacity.

Our range of optional equipment means you can customise your filter to meet your requirements.



#### BLASTONE INTERNATIONAL (AUST) PTY LTD

57 Alexandra Place Murarrie Queensland 4172 Australia

PO Box 3385 Tingalpa DC Queensland 4172 Australia

Toll Free 1800 190 190 Tel (07) 3329 5000 Fax (07) 3329 5066 Email sales.au@blastone.com www.BlastOne.com

#### AUSTRALIAN OFFICES

Adelaide | Brisbane | Darwin | Mackay | Melbourne | Perth | Sydney

#### BLASTONE INTERNATIONAL

4510 Bridgeway Avenue Columbus Ohio 43219 USA

Toll Free 800 999 1881 Tel (614) 476 3000 Fax (614) 476 3002 Email sales@blastone.com www.BlastOne.com

#### NORTH AMERICAN OFFICES

Chicago | Columbus | Los Angeles | Minneapolis